

# B-100 Series Ultraviolet Lamps

## Operating Instructions

---


These Operating Instructions cover UVP's high-intensity UV lamp models (115V and 220V):

**B-100AP • B-100AP/R • B-100YP • B-100Y • B-100A • B-100A/R**

The B-100AP, B-100AP/R and B-100YP models are manufactured with UVP's Cool-Touch™ plastic housing. B-100A, B-100A/R and B-100Y models are manufactured of a durable scratch-resistant, powder-painted aluminum. The B-100Y and B-100YP models feature a yellow filter. Lamp and accessory drawings and replacement parts lists for each model are included in these instructions.

### General Information for Ultraviolet Inspection

Make tests in a darkened area for best contrast. Fluorescence may be bright or dim, depending on amount of fluorescing agent present. Be careful to distinguish fluorescence from blue light reflection on polished metal surfaces. UVP Model UVC-503 Safety Goggles are recommended to eliminate blue haze resulting from eye fluorescence caused by reflected longwave ultraviolet radiation.

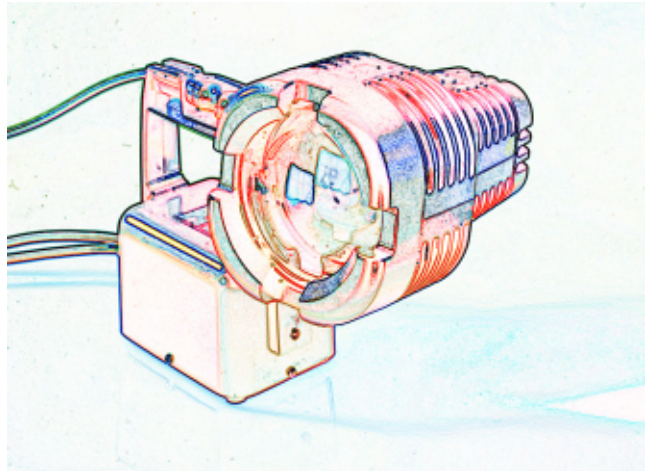

### Installation

An 8-foot primary and secondary line cord is supplied on all units except the B-100AP/R and B-100A/R which come with 8-foot primary and 20-foot secondary line cords. If the receptacle is not fitted with 3-prong plug, use a 2-prong adapter and attach grounding lead to screw center of receptacle plate.

### Operation

Turn on/off switch on transformer base to ON position for starting. Although a cool lamp starts immediately, allow approximately five minutes to reach full brilliance. When use is intermittent, it is preferable to leave the lamp burning. This assures longer bulb life and avoids delay in restarting. Restarting does not take place immediately as the mercury vapor within the bulb must cool. Delay is typically 4 to 8 minutes. Any interruption in the power supply, such as an excessive drop in voltage, will extinguish the arc. Restarting typically will be delayed while the mercury vapor cools. Full voltage must be restored to permit the arc to strike.

While mounted on the transformer base, the lamp housing can rotate 360°. For handheld use, detach the lamp housing from the transformer.


## Replacement Parts List

Parts lists are provided for the B-100 lamp models. Two parts lists are shown -- one for the Cool-Touch housing lamps and one for the aluminum housing lamps. For additional ordering information or an authorized dealer, contact UVP, Inc. in Upland, California at (800)452-6788 or (909)946-3197; Ultra-Violet Products Limited in England at +44(0)1223 420022.

### B-100AP, B-100AP/R AND B-100YP LAMPS


115 VOLT, 60 Hz; 220 VOLT, 50 Hz

PART NO.	PART NAME
34-0054-01	SPOT BULB
36-0043-02	SOCKET, MEDIUM BASE
38-0168-01	BEZEL AND UV FILTER ASSEMBLY
38-0168-02	BEZEL & YELLOW FILTER ASSEM.
88-0016-01	BALLAST REPLACEMENT KIT, 115V
88-0016-02	BALLAST REPLACEMENT KIT, 220V
77-0007-01	B-100AP LAMP HEAD ASSEMBLY, 115V
77-0007-04	B-100AP LAMP HEAD ASSEMBLY, 220V
77-0007-05	B-100AP LAMP HEAD ASSEMBLY, 115V
77-0007-06	B-100AP/R LAMP HEAD ASSEMBLY, 220V
77-0007-07	B-100YP LAMP HEAD ASSEMBLY, 115V
77-0007-08	B-100YP LAMP HEAD ASSEMBLY, 220V

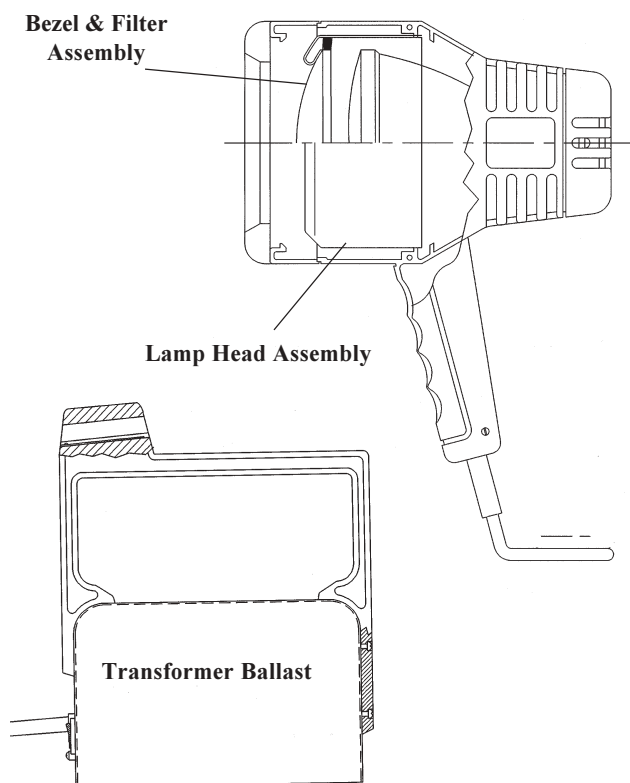


## Accessories

The **Lamp Funnel**, for use with the Cool-Touch models, can be used to direct the UV source directly where you want it, eliminating stray UV radiation. To attach, remove plastic bezel from unit by pulling tabs apart and sliding bezel out. Next insert retaining clips in the four slots of the interior bezel as shown in Section A-A. Hold funnel in place and reattach assembled bezel, over funnel, to unit by pulling tabs apart and snapping into place.

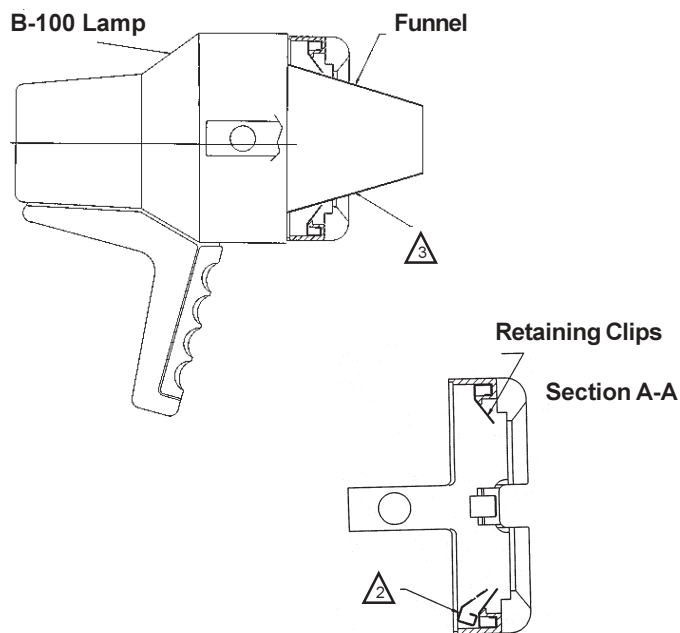

PART NO.	PART NAME
19-0126-01	LAMP FUNNEL (for B-100AP, YP models)

The **Finger Guard**, to keep fingers away from the filter glass, easily snaps into place on the face of the Cool-Touch lamps.

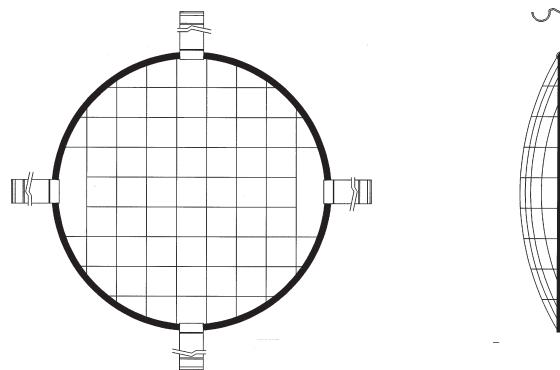

PART NO.	PART NAME
89-0068-01	LAMP FUNNEL (for B-100AP)


## B-100A, B-100A/R AND B-100Y LAMPS

115 VOLT, 60 Hz; 220 VOLT, 50 Hz

(Note: B-100A/R NSN 6635-611-5617 ordering number)

PART NO.	PART NAME
34-0054-01	SPOT BULB
36-0043-02	SOCKET, MEDIUM BASE
38-0168-03	BEZEL & UV FILTER ASSEM.
38-0168-04	BEZEL & YELLOW FILTER ASSEM.
77-0005-01	LAMP HEAD ASSEMBLY WITH 8ft. CORD, 115V
77-0005-02	LAMP HEAD ASSEMBLY WITH 8ft. CORD, 220V
77-0005-03	LAMP HEAD ASSEM. WITH 20ft. CORD, 115V
77-0005-04	LAMP HEAD ASSEM. WITH 20ft. CORD, 220V
88-0016-01	BALLAST REPLACEMENT KIT, 115V
88-0016-02	BALLAST REPLACEMENT KIT, 220V
19-0124-01	LAMP VISOR



The **Lamp Visor** shields the user's eyes from the UV and directs the UV on the object of analysis. Visor shown below is used with Cool-Touch lamps. Visor clips onto the top of the lamp head. Front and side facing views shown.

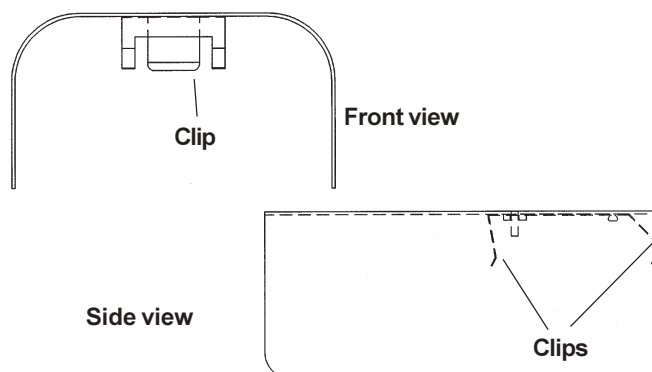

PART NO.	PART NAME
19-0125-01	LAMP VISOR (for B-100AP, YP models)

## Accessories

The **Lamp Visor** shields the user's eyes from the UV and directs the UV on the object of analysis. Visor clips onto the sides of the lamp head. Side and front facing views shown.

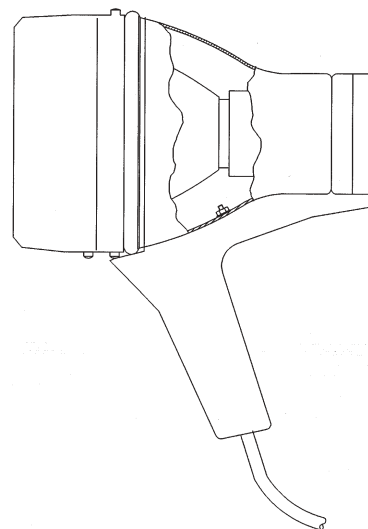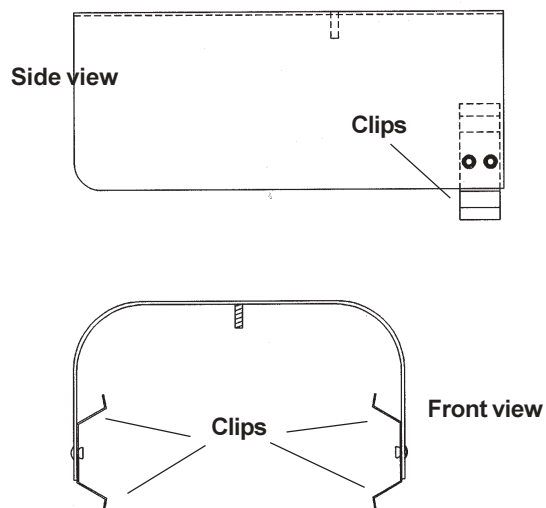

## Bulb Replacement and Repairs

- To replace the bulb, first remove the filter assembly:
  - Expand *front bezel retainers* over circular plastic holders and slide front bezel forward and off.
  - Expand the two silver filter assembly spring clips outward to release filter assembly from lamp body.
- The bulb can be replaced after filter assembly has been detached. When replacing the bulb, screw bulb in tightly to prevent loosening from heat expansion and contraction during operation.

Prior to returning any lamps to the UVP factory for repair, contact UVP's customer service department to obtain a **Return Goods Authorization (RGA)** number.

## Bulb Standardization

For testing uniformity, a standard bulb intensity is required. A special light meter calibrated for measuring ultraviolet intensity is required to check bulb standardization.

The UVX Radiometer and UVX-36 sensor or the J-221 UV Meter employ a dense ultraviolet filter and is specifically calibrated to read ultraviolet intensity of 365nm in mW/cm<sup>2</sup> from a line source lamp. An ordinary light meter will not deliver readings of ultraviolet energy.

## CAUTION

1. Wear hand and eye protection if filter glass is cracked. Replace broken filter as soon as possible.
2. To avoid electrical shock hazard, unplug lamp before servicing.
3. Bulb runs very hot; allow unit to fully cool before servicing.
4. Bulb contains mercury. If bulb breaks, wear gloves to remove.
5. B-100A, B-100A/R and B-100Y lamp housings become hot during operation and care should be taken not to touch them. Lamp handles do not heat even after several hours of use.

## PRODUCT WARRANTY

UVP's products are guaranteed to be free of defects in materials, workmanship and manufacture for one (1) years from the date of purchase; transilluminators are guaranteed for (2) years. Consumable and disposable parts including, but not limited to bottles, tubes and filters, are guaranteed to be free from defects in manufacture and materials for ninety (90) days from date of purchase. If equipment failure or malfunction occurs during the warranty period, UVP shall examine the inoperative equipment and have the option of repairing or replacing any part(s) which, in the judgment of UVP, were originally defective or became so under conditions of normal usage and service.

No warranty shall apply to any instrument, or part thereof, that has been subject to accident, negligence, alteration, abuse or misuse by the end-user. Moreover, UVP makes no warranties whatsoever with respect to parts not supplied by UVP or that have been installed, used and/or serviced other than in strict compliance with the instructions appearing in the operational manual supplied to the end-user.

In no event shall UVP be responsible to the end-user for any incidental or consequential damages, whether foreseeable or not, including, but not limited to property damage, inability to use equipment, lost business, lost profits, or inconvenience arising out of or connected with the use of instruments produced by UVP. Nor is UVP liable or responsible for any personal injuries occurring as a result of the use, installation and/or servicing of equipment. This warranty does not supersede any statutory rights that may be available in certain countries.



ULTRA-VIOLET PRODUCTS

**Corporate Headquarters: UVP, Inc.** 2066 W. 11th Street, Upland, California 91786  
Tel: (909)946-3197 or Toll Free (800)452-6788 • Fax: (909)946-3597 • E-mail: [uvp@uvp.com](mailto:uvp@uvp.com)

**European Operations: Ultra-Violet Products Ltd,** Unit 1, Trinity Hall Farm Estate, Nuffield Road, Cambridge, CB4 1TG UK  
Tel: +44(0)1223 420022 • Fax: +44(0)1223 420561 • E-mail: [uvp@uvp.co.uk](mailto:uvp@uvp.co.uk)